

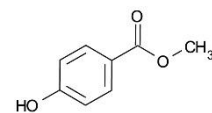


## Conditions


Column: ACE 3 Phenyl  
Dimensions: 150 x 2.1 mm  
Part Number: ACE-115-1502  
Mobile Phase: 25 mM ammonium acetate pH 6.8 in H<sub>2</sub>O/MeOH (50:50 v/v)  
Flow Rate: 0.2 mL/min  
Injection: 2 µL  
Temperature: 40 °C  
Detection: UV, 240 nm

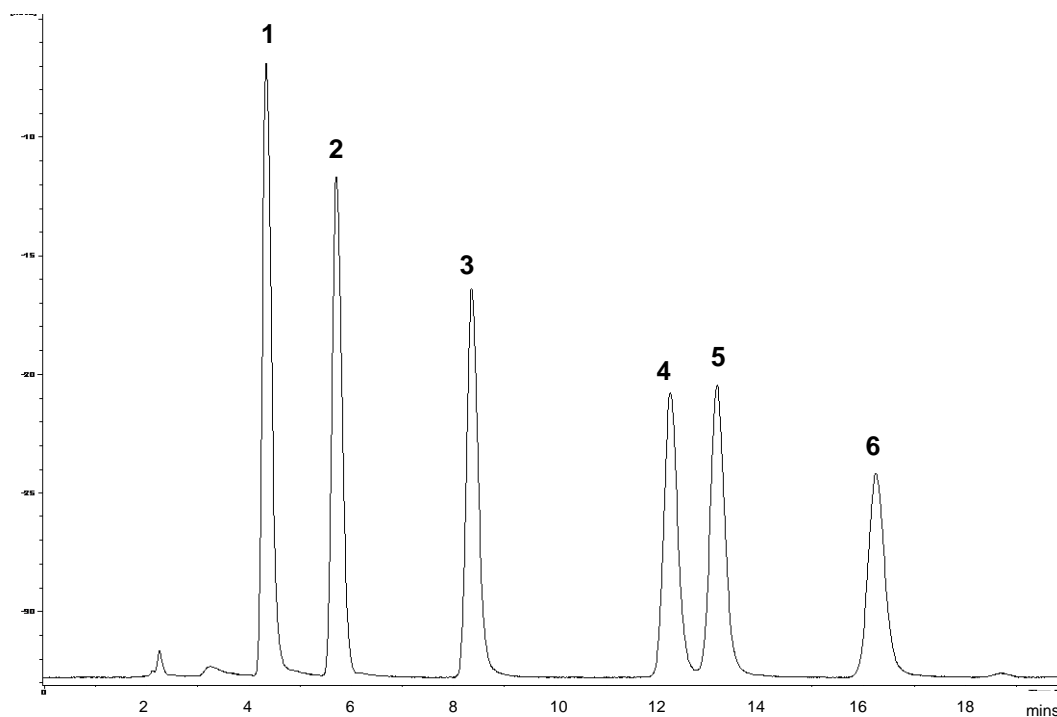


1. Methyl paraben

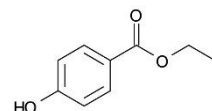


2. Ethyl paraben

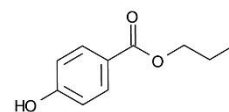


3. n-Propyl paraben

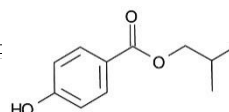

4. i-Butyl paraben

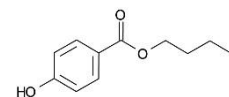

5. n-Butyl paraben

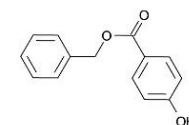

6. Benzyl paraben