

# Rapid analysis of Cold and Flu Medicine Components

**ACE**<sup>®</sup>  
Ultra-inert  
UHPLC & HPLC Columns


Application #AN4780

## Conditions

Column: ACE Excel 2 C18-AR  
Dimensions: 100 x 3.0 mm  
Part Number: EXL-109-1003U  
Mobile Phase: A: 0.1% formic acid in H<sub>2</sub>O  
B: 0.1% formic acid in MeCN

Time (mins)	%B
0	5
3.6	20
5.4	80
6	80
6.5	5

Flow Rate: 1.0 mL/min  
Injection: 5 µL  
Temperature: 40 °C  
Detection: UV, 275 nm  
Instrument: Chromaster UltraRs

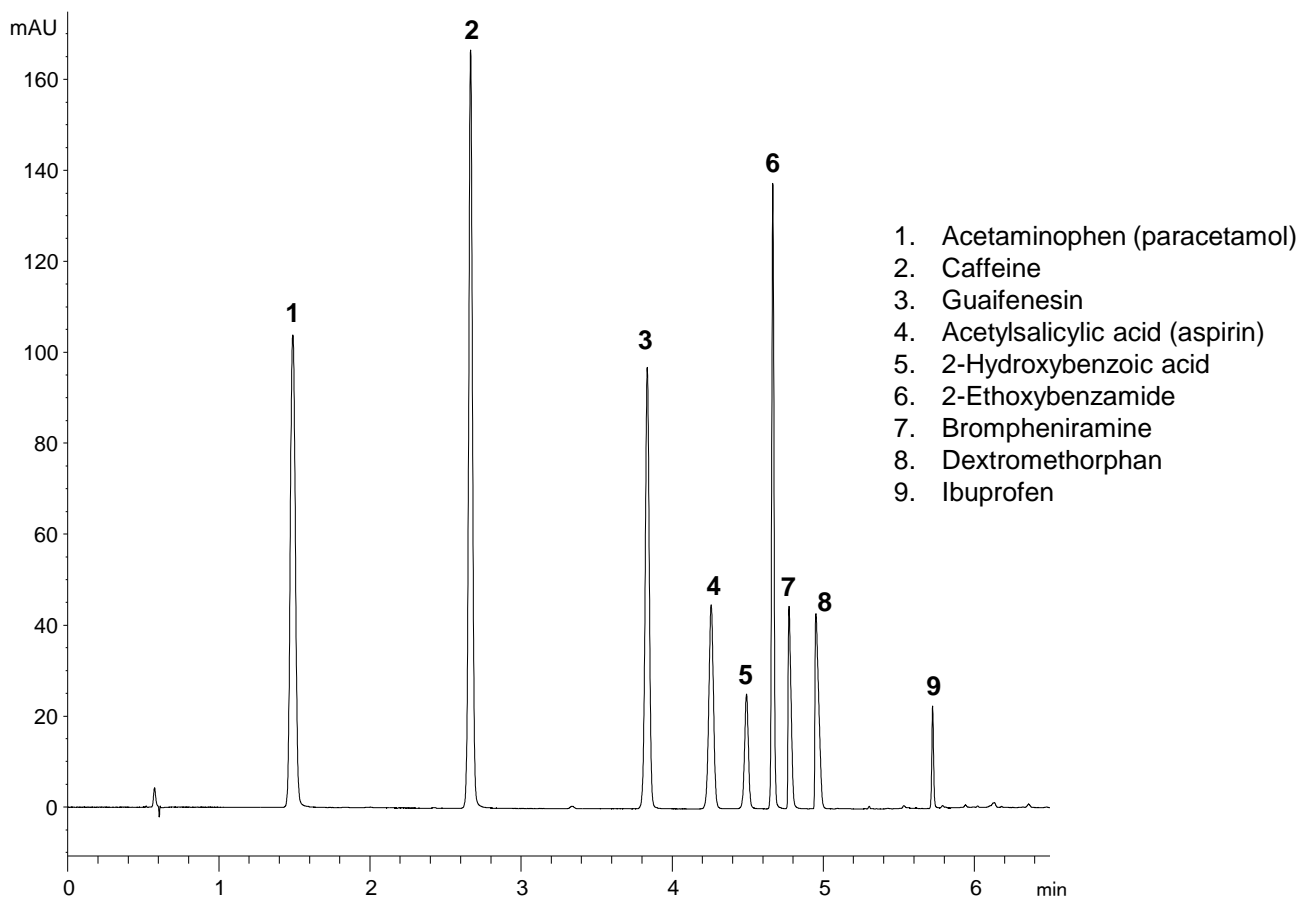

© 2017 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit  
[www.ace-hplc.com](http://www.ace-hplc.com) or email: [info@ace-hplc.com](mailto:info@ace-hplc.com)

  
**ACE**<sup>®</sup>  
UHPLC & HPLC  
Columns